

CLINITEX[®] HealWell[®]

Night Splint

with  **Microban[®]**
antimicrobial protection

Model No. 58-409

Easy Push Lock Buckles

Removable Liner -

Is made of the unique C3™ Tri-Laminate Material

Adjustable Bilateral Dorsiflexion Straps -


Easy to adjust tension straps apply a controlled amount of Dorsiflexion (0 -10°) and stretch to the Plantar Fascia

Removable Foam Wedge -

Provides additional toe extension

Non-Slip Bottom -

Not intended for Ambulation but helps prevent accidental slipping on flooring

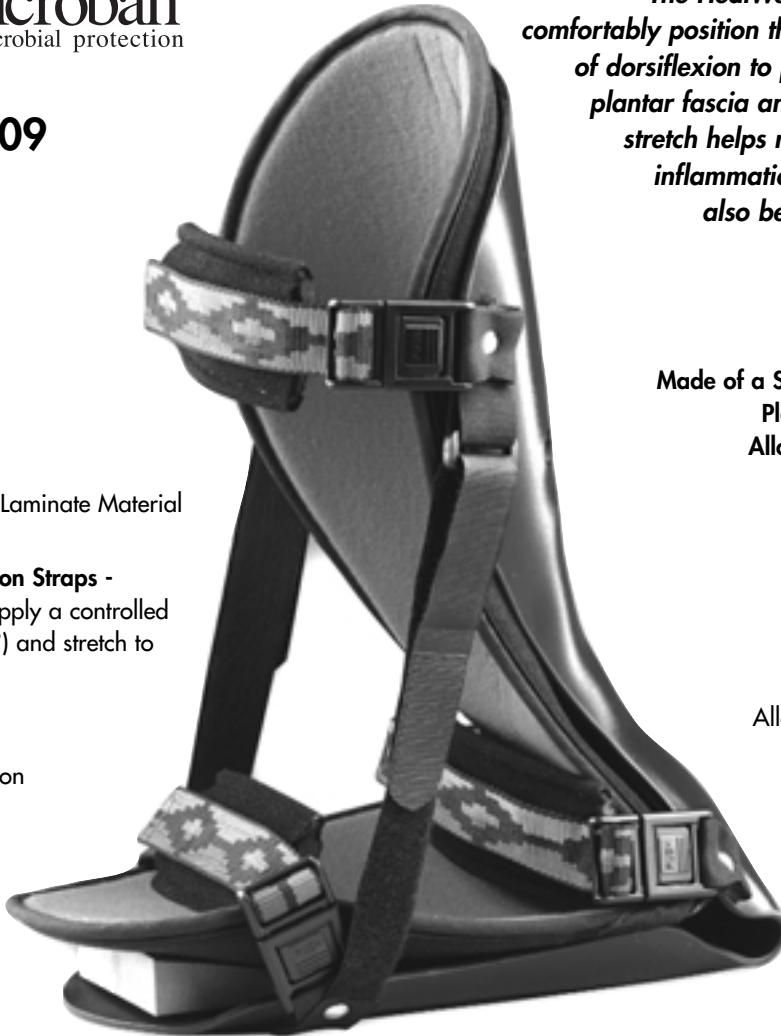

Indicated to Treat:
Plantar Fasciitis • Achilles Tendinitis
Plantarflexion Contractures • Heel Spurs
Equinus Conditions • Pronation Syndromes
Muscle Tightening and Runner's Cramp

*The HealWell[®] Night Splint is designed to comfortably position the foot in a controlled amount of dorsiflexion to provide a gentle stretch to the plantar fascia and Achilles Tendon. This gentle stretch helps reduce the muscle contracture, inflammation and associated pain. It has also been clinically proven to reduce symptoms of plantar fasciitis.**

Made of a Spring-Action Lightweight Kydex[®] Plastic with Ventilated Calf Area to Allow for Comfort and Breathability

Three Strategically Placed Comfortable Padded Straps -
Ensure the foot is comfortably maintained in the correct position

Open Cut Heel Design -
Allows for the foot to be comfortably and properly positioned

Strong Kydex[®] Plastic -
Prevents unwanted torquing of splint at the ankle

Patent No. D434,504



SIZING CHART


Size	Men's Shoe Size	Women's Shoe Size
Small	Under 6	5 - 7
Medium	6 - 8	7 1/2 - 9 1/2
Large	8 1/2 - 10 1/2	10 - 12
X-Large	11 - 13	Over 12
XX-Large	Over 13	

Discover the  **Tri-Laminate** benefits for easy and comfortable wearing.

OUTER LAYER: Soft, hook compatible outer laminate is extremely durable and comfortable. Black color hides dirt for easy and attractive long term wear.

MIDDLE LAYER: Foam inner is treated with Microban[®] Antimicrobial Protection. The Microban[®] agent is actually part of the foam's molecular structure and helps provide continuous control over the growth of certain odor and stain causing bacteria, yeasts and fungi. **The result: A Clean and Comfortable product.**

INNER LAYER: Soft liner is made of a unique material that transports moisture away from the skin and into the Microban[®] treated foam layer. **The result: Dry, Cool and Comfortable wear.**

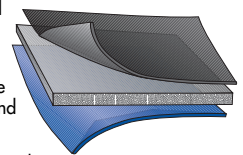

Suggested HCPCS Code: L4396

©2000 FLA Orthopedics, Inc.

FLA ORTHOPEDICS[®]

1-800-327-4110

www.HealWellinfo.com

* "The Use of Night Splints for Treatment of Recalcitrant Plantar Fasciitis." Foot & Ankle, Vol. 12 No. 3, Dec. 1991., Keith L. Wapner, M.D. and Peter Sharkey, M.D. Philadelphia, Pennsylvania and Oakland, California

Microban[®] is the Registered Trademark of the Microban Products Company, Huntersville, NC
C3[™] is the Trademark of FLA Orthopedics, Inc., Miramar, FL
HealWell[®] is the Registered Trademark of FLA Orthopedics, Inc., Miramar, FL